

Field Service Procedure

Fabius Tiro™ Pressure Control Upgrade Anesthesia



Revision -
1/21/04
SP00302

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

Service Procedure

Fabius Tiro Pressure Control Upgrade

FABIUS TIRO PRES- SURE CONTROL UPGRADE

1 SWITCH PANEL INSERT REPLACE- MENT

**NOTE:**

The screen illustrations contained in this document are for reference only and therefore may or may not reflect the current equipment, software version being installed or exact Fabius product.

1. Set the main power switch to OFF and disconnect AC power from the machine.
2. Swing open the ventilator door.
3. While supporting the bezel assembly, remove the two screws securing the left side of the bezel.

**CAUTION:**

Do not pull the bezel outward as this will break the locking tabs at the right side of the bezel.

4. Slide the bezel to the right to release its locking tabs from the machine. See [Fig. 1](#).

**NOTE:**

Verify that the bezel locking tabs are clear of the machine before pulling the bezel away from the machine.

Fabius Tiro Pressure Control Upgrade

Service Procedure

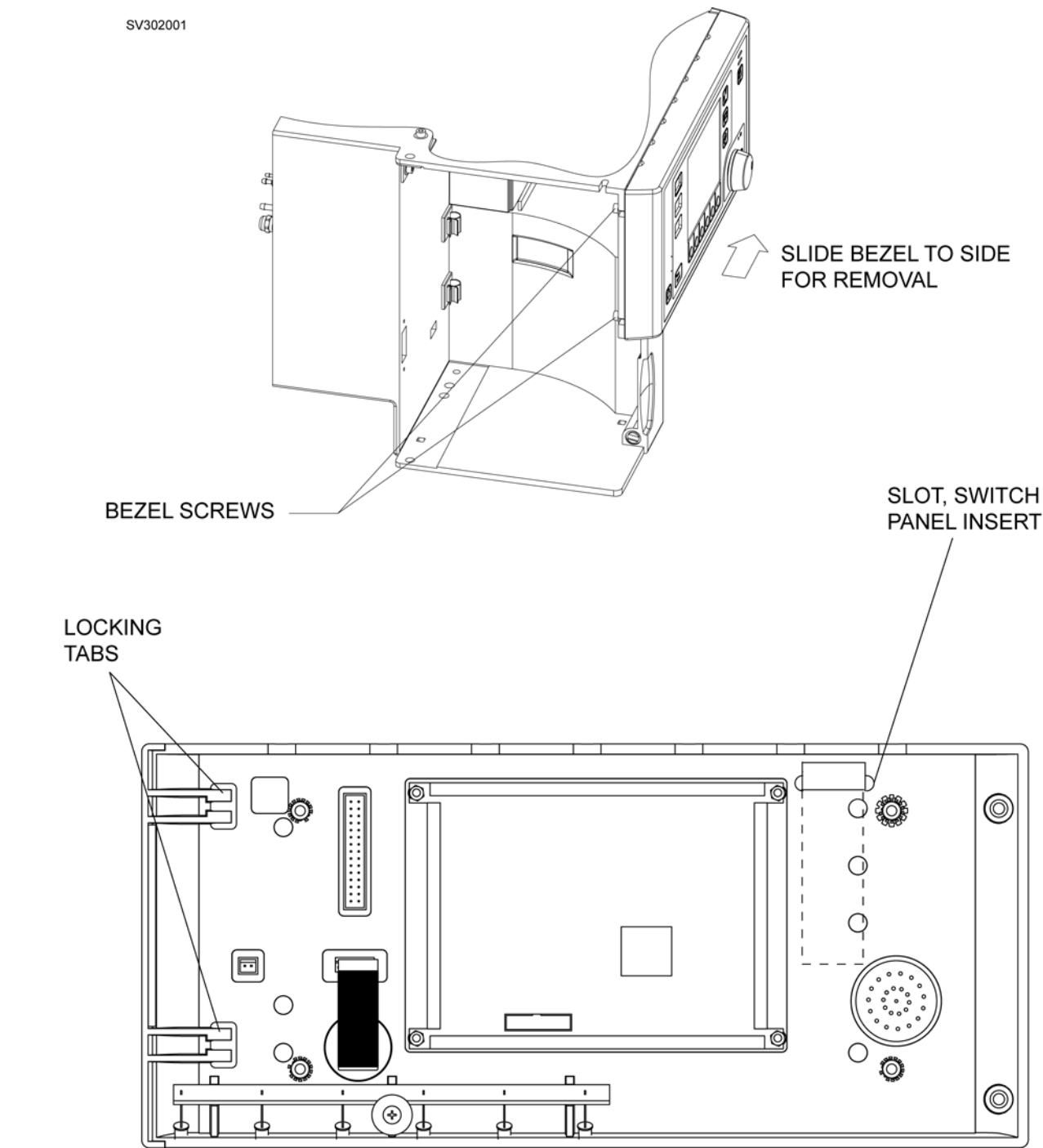


Fig. 1 Bezel Removal and Switch Panel Insert Slot Location

Service Procedure

Fabius Tiro Pressure Control Upgrade

5. At the Switch Panel Insert Slot, locate the exposed end of the current switch panel insert and pull to remove it from the bezel assembly. See [Fig. 1](#).
6. Select the correct switch panel insert per local application or customer requirement. Refer to [Table 1](#) for corresponding part numbers.

Table 1 Switch Panel Insert Identification

PART NUMBER	ALTERNATE PART NUMBER	DESCRIPTION	QTY
4118120-002	4199860	Insert, SW Panel - INTL V/P	1
4118120-006	4199856	Insert, SW Panel - DOM V/P	1
4118120-010	4199852	Insert, SW Panel - FR V/P	1

7. With the replacement switch panel insert correctly oriented, insert into the bezel assembly at the switch panel insert slot. See [Fig. 1](#).
8. Verify the replacement switch panel insert lettering is centered and level within each respective keypad.
9. Carefully fold the cables behind the bezel assembly as you bring it into position on the machine. Hold the bezel against the machine and slide it to the left to engage the locking tabs in the slots in the machine housing.
10. While supporting the bezel assembly, reinstall the two screws at the back of the bezel on the left side.
11. Close the ventilator door.
12. Restore power to the machine and verify all of the keypad and display functions.
13. Proceed with Secure Options Section.

2 SECURE OPTIONS



NOTE:

Secure Options screens are located in a sub-menu of the main Service Menu on the host Fabius Tiro units with software version 2.N.
PATH: Main Service Screen/Configure/Secure Options

1. Access the Secure Options screen located in the Main Service Screen and obtain the Device ID#, and (Option) Description and Part Number. For Pressure Control, see "PCV Vent Mode".
2. Contact DrägerService - Technical Support: Phone 1-800-4-Dräger, Phone: 215-721-5402, or e-mail to techsupport@draegermed.com and provide the following information to receive the necessary release code(s):
 - Machine Type (Fabius GS or Fabius Tiro)
 - (Option) Description and Part Number
 - Machine Serial Number
 - Device ID Number

Fabius Tiro Pressure Control Upgrade

Service Procedure

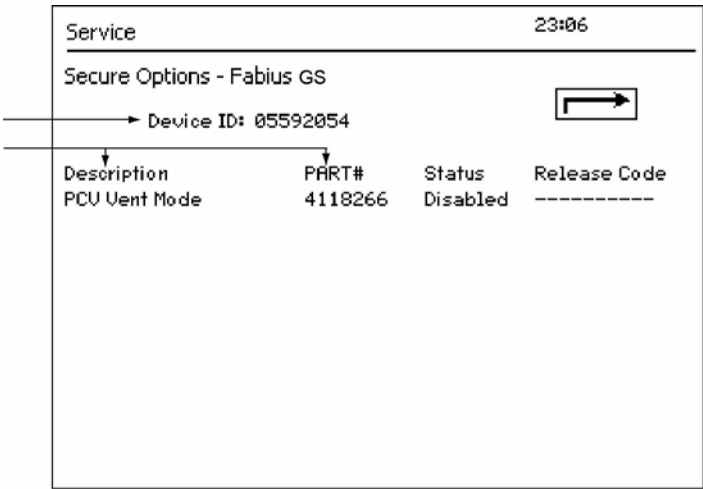


Fig. 2 Secure Options Screen

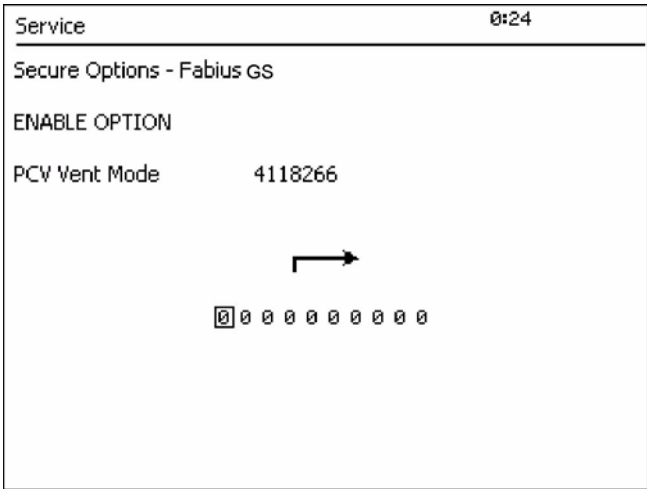


Fig. 3 Secure Options Screen

3. Access the Secure Options screen and enter the release code number provided by Technical Support and confirm. This enables the ventilator option(s).
4. From the Standby screen, press the 'Pressure Control' key and confirm. Verify pressure control is active, then press 'Standby' key and confirm.
5. Perform all PMS test steps contained in Section 8.20 (Ventilator) of the Fabius Tiro Service Manual to verify ventilator functions.

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

In the United States:
Draeger Medical, Inc.
3136 Quarry Road
Telford, PA 18969
USA

Tel: (215) 721-5402
(800) 543-5047
Fax: (215) 721-5785

Fabius Tiro Pressure Control Upgrade en
SP00302

Subject to change without notice.
Will not be replaced in the event of modifications.

In all other countries:
Dräger Medical AG & Co. KGaA
Moislinger Allee 53 – 55
D-23542 Lübeck
Germany

Tel: (++49) (0) 1805-3723437
Fax: (++49) 451/882 - 3779